



COMBAT MEDICAL PRODUCTS



GeM
Government
e Marketplace



"Saving life of a soldier is more important than killing an enemy"



FIRST FIELD DRESSING

As per DGQA specification no.
Prov/S/4665/TC-1(a)/2016/01/(a)

Recommended for:

To provide compression
on bleeding point of the limbs

The **EmRescue FFD** is a specialized compression bandage designed for quick deployment by soldiers in the event of a severe injury where halting blood loss is of utmost importance. Its distinctive design allows it to be self-applied by an injured soldier using just one hand. The FFD generates adequate pressure to halt bleeding and can maintain that pressure on the wound until the injured soldier receives comprehensive medical care.

Key Features:

- It includes a medical-grade pressure-sensitive adhesive at the beginning of the dressing, facilitating single-handed bandage application.
- The wound pad is constructed with non-leachable silver, which is bactericidal, or it can be impregnated with Chlorohexidine Di Gluconate as an antimicrobial treatment. This special wound pad not only prevents bacterial growth but also actively eliminates bacteria, reducing the risk of wound infection.
- The wound pad is comprised of a super-absorbent material capable of absorbing fluid up to nearly 50 times its weight. This high absorption capacity allows the dressing to effectively absorb and retain a significant quantity of blood.
- To prevent accidental unfolding, holding knots are incorporated between the layers of the bandage. This ensures that the bandage remains securely in place even if dropped on the ground.
- A specially designed ridge is incorporated into the FFD, which applies pressure to the bleeding point. This design reduces the risk of gangrene, as blood circulation can be maintained around the sides of the ridge, which remain pressure-free.
- Due to its high elasticity, the FFD can accommodate post-injury swelling, ensuring a comfortable fit.
- It includes a hook and loop fastener (velcro) at the end of the dressing, making application easy and eliminating the need for additional materials to secure the bandage.
- To uphold sterility, the FFD is double-packed and supplied in a sterilized condition. The inner packaging is made of medical-grade paper and polyester laminate, while the outer packaging consists of transparent polyester laminate sheet.
- The FFD is vacuum-packed in a robust aluminum packaging, offering protection against light, ultraviolet radiation, water vapor, oxygen, and microorganisms, ensuring its integrity during storage and use.

Instructions for use



1 Teas open the pouch and take out the inner pack. Peel open inner pack & take out the dressing.



2 Remove release liner from end.



3 Holding knots have been provided between layers of bandage so that bandage would not unfold on its own even when it is dropped on ground



4 Press the end before the wound. The ridge must oppose the bleeding wound. Wrap bandage around the limb.



5 Press the free flap '3B' against '3A'



6 Pull the bandage forcefully and wrap it around the limb maintaining strong pull




7 Keep on wrapping the bandage tightly.



8 End part of bandage is cohesive. Press it on the bandage and it will stick on it.

Ordering information	First Field Dressing
Ordnance Item Part No.	KND NIV CL2 248 FFD
Product Code	FFD
NSN No.	6510720499290
GeM Catalogue ID	5116877-57747712587-cat
GeM Product ID	7751183-15165600236

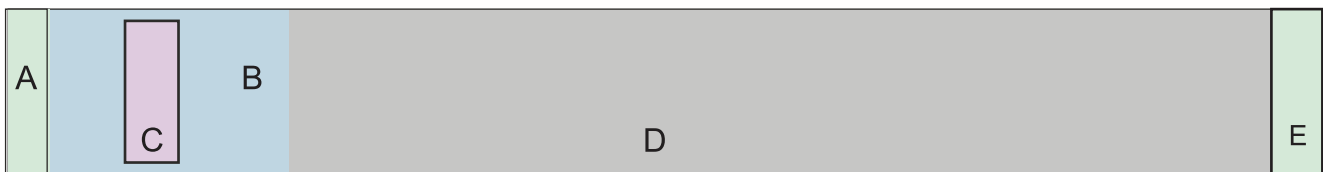

A: Adhesive Strip B: Super Absorbent Wound Pad with Non Leachable Silver/ Anti Septic Wound Pad
 C: Rigid Ridge D: Main Bandage E: Velcro Strip



Made in India by:

G. Surgiwear Limited
Rasoolpur, Jahanganj, Shahjahanpur-242001 (INDIA)
Tel: +91 7518700177 Fax: +91 7518700170
e.mail: info@surgiwear.in
CIN : U24236UP1990PLC012073

Exclusive Distributor:

Medicaid Agencies
UG - 4&5, Shivlok House II
Karampura Commercial Complex
New Delhi - 110015 (INDIA)
Phone: +91-11-41427278, 25920509
Mob. : +91- 9810023232
e.mail : info@medicaidagencies.com

Manufacturing Quality Products for

Neurosurgery	Cardiothoracic
Orthopedics	Laparoscopic
Ophthalmology	Obst. Gynae
Plastic Surgery	Burn Treatment
Urosurgery	Wound Dressing
Andrology	Emergency
General Use	Pressure Dressing

* G Surgiwear Limited retains the authority to make adjustments, modifications, additions, or removals to the products listed in the brochure as deemed necessary, without prior notice. Product specifications are subject to change without prior notice. **EmRescue is registered brand of G Surgiwear Limited.